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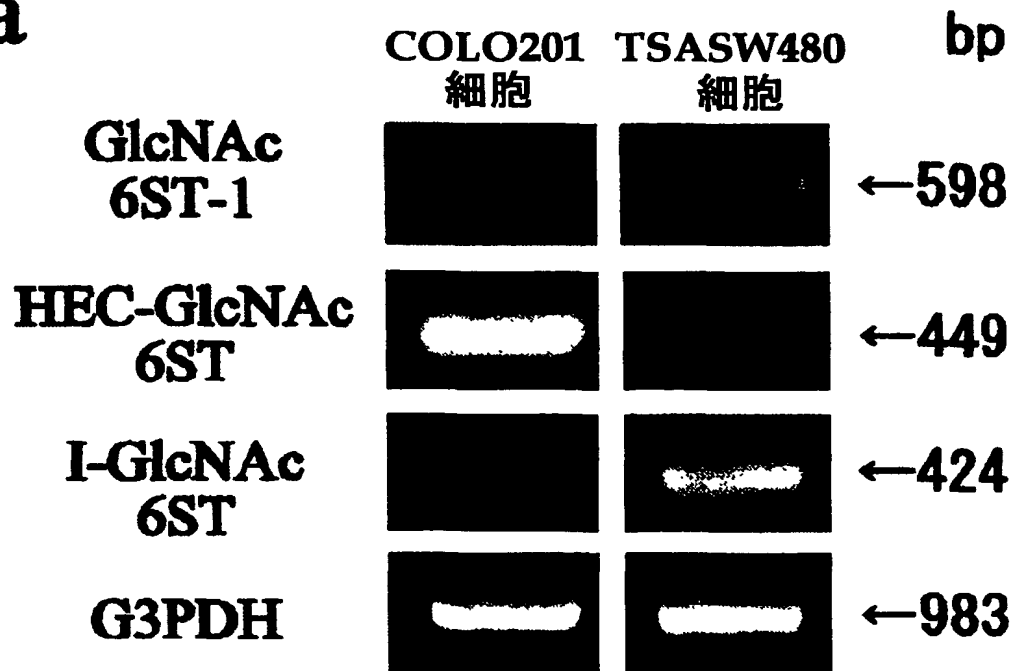
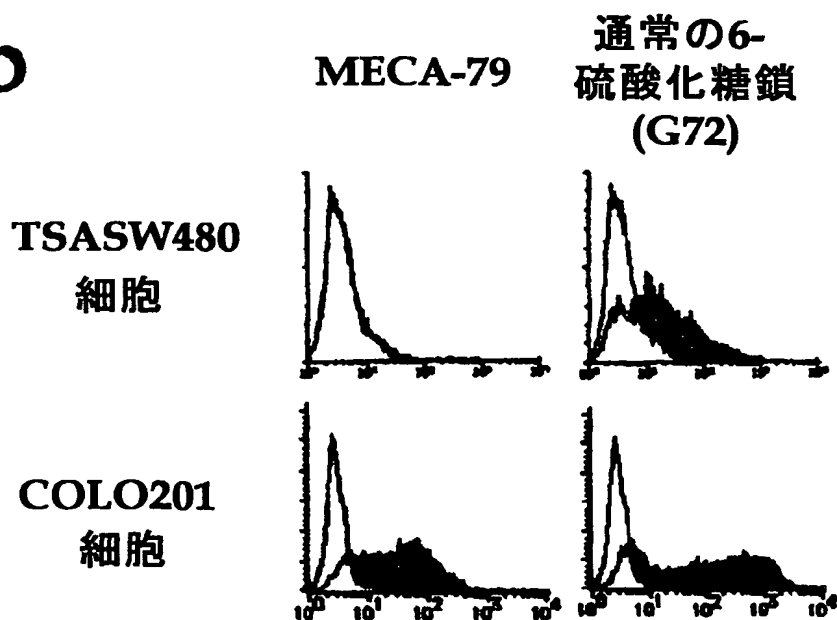
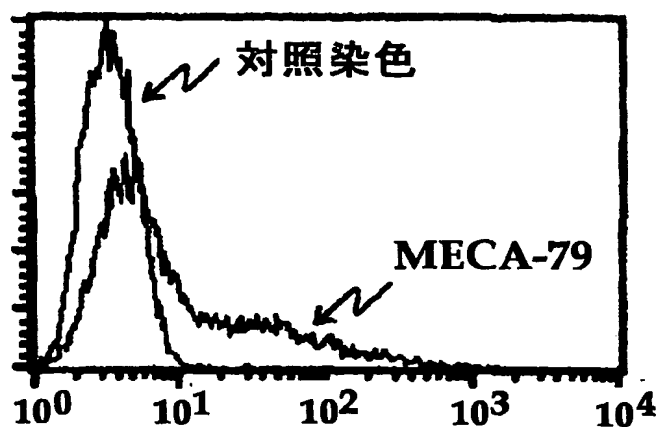
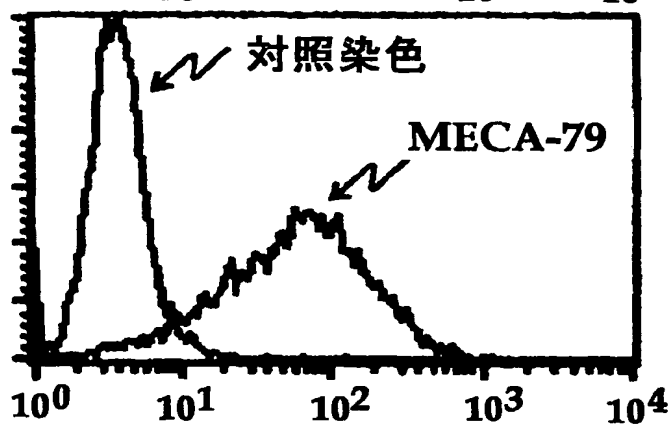
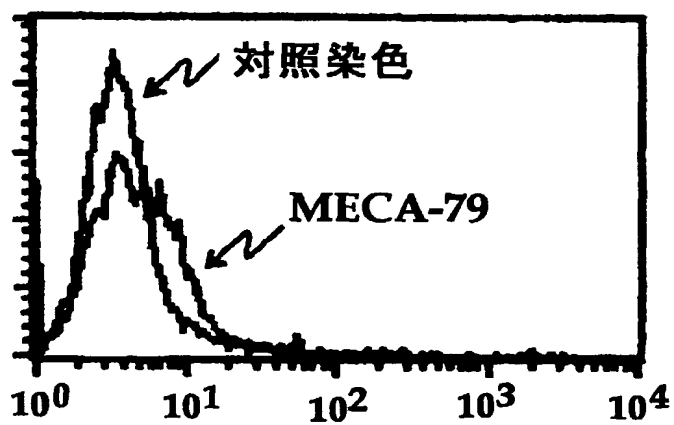
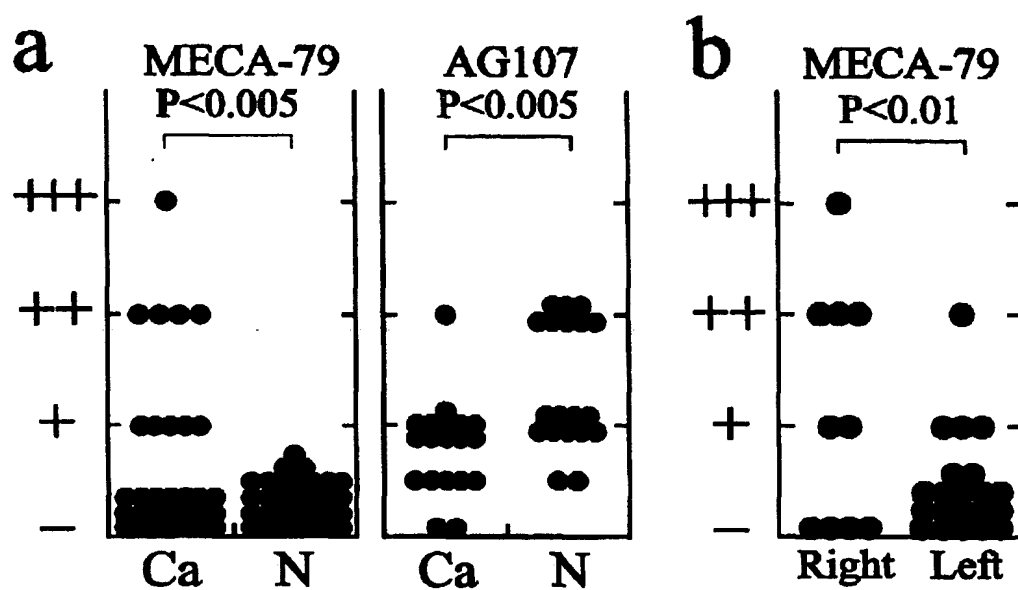
BD PharMingen Technical Data Sheet

Page 1 of 2

PURIFIED RAT ANTI-MOUSE PNAd CARBOHYDRATE EPITOPE (CD62L Ligand) MONOCLONAL ANTIBODY**PRODUCT INFORMATION**

Catalog Number: **553863 (Was: 09961D), 0.5 mg**
 Description: Purified anti-mouse PNAd Carbohydrate Epitope (CD62L Ligand)
 Clone: MECA-79
 Immunogen: Collagenase-dispersed BALB/c lymph node stroma¹
 Isotype: Rat (Wistar) IgM, κ
 Contents: Purified immunoglobulin in 10 mM phosphate buffer, pH 7.2 with 500 mM NaCl and 0.09% (w/v) sodium azide.

SPECIFICITY

The MECA-79 antibody reacts with sulfate-dependent carbohydrate epitopes of peripheral lymph node addressin (PNAd).² The MECA-79-reactive antigen is closely associated with the carbohydrate ligands for L-selectin (e.g., CD34, GlyCAM-1, MAdCAM-1), which are expressed on high endothelial venules (HEV) in lymphoid tissues and at sites of chronic inflammation.^{1,2,3,4,5,6} Cross-reactivity with human,^{3,4} ovine,⁷ bovine,⁷ primate,⁷ and porcine⁸ tissues has been observed. MECA-79 antibody inhibits L-selectin-dependent lymphocyte and platelet homing to lymph nodes *in vivo*.⁹ and *in vitro* adhesion to lymphoid tissue HEV^{1,4} and immobilized PNAd.^{3,9,10}

PREPARATION AND STORAGE

The antibody was purified from tissue culture supernatant by affinity chromatography. The antibody solution should be stored undiluted at 4°C.

USAGE

This antibody has been tested by immunohistochemical staining (IHC) of citrate-pretreated formalin-fixed paraffin-embedded sections (5 - 20 μ g/ml) to assure specificity and reactivity. Other reported applications include IHC of acetone-fixed frozen sections,^{1,4,5} immunoprecipitation,^{2,3} western blot analysis,¹⁰ and *in vitro* and *in vivo* adhesion blocking.^{1,3,4,9,10} Since applications vary, each investigator must determine dilutions appropriate for individual use.

Caution: Sodium azide is a reversible inhibitor of oxidative metabolism; therefore, antibody preparations containing this preservative agent must not be used in cell cultures nor injected into animals. Sodium azide may be removed by washing stained cells or plate-bound antibody or dialyzing soluble antibody in sodium azide-free buffer. Since endotoxin may also affect the results of functional studies, we recommend the NA/LE™ (No Azide/Low Endotoxin) antibody format for *in vitro* and *in vivo* use.

REFERENCES

1. Streeter, P.R., B.T.N. Rouse, and E.C. Butcher. 1988. Immunohistologic and functional characterization of a vascular addressin involved in lymphocyte homing into peripheral lymph nodes. *J. Cell Biol.* 107: 1853 - 1862.
2. Hemmerich, S., E.C. Butcher, and S.D. Rosen. 1994. Sulfation-dependent recognition of high endothelial venules (HEV)-ligands by L-selectin and MECA 79, an adhesion-blocking monoclonal antibody. *J. Exp. Med.* 180: 2219 - 2226.
3. Berg, E.L., M.K. Robinson, R.A. Warnock, and E.C. Butcher. 1991. The human peripheral lymph node vascular addressin is a ligand for LECAM-1, the peripheral lymph node homing receptor. *J. Cell Biol.* 114: 343 - 349.
4. Michie, S.A., P.R. Streeter, P.A. Bolt, E.C. Butcher, and L.J. Picker. 1993. The human peripheral lymph node vascular addressin. An inducible endothelial antigen involved in lymphocyte homing. *Am. J. Pathol.* 143: 1688 - 1698.
5. Faveeuw, C., M.-C. Gagnerault, and F. Lepault. 1994. Expression of homing and adhesion molecules in infiltrated islets of Langerhans and salivary glands of nonobese diabetic mice. *J. Immunol.* 152: 5969 - 5978.

Please see Page 2.

BD Biosciences www.bdbiosciences.com**Biometric Imaging
Clontech**

Alexis Pacific
Tel (530) 861-6323
Fax (530) 860-1590

Brazil
Tel (55) 11 5102 9831
Fax (55) 11 5105 9834

Canada
Tel (800) 599-0187
Fax (905) 543-8391
canada@bdb.com

**Immunocytometry Systems
Labware**

Europe
Tel (33) 33 720 111
Fax (33) 33 720 450

Europe
Tel (408) 6321 805 521
Fax (408) 6321 305 531

**Transduction Laboratories
PharMingen**

Japan
Tel (81) 3 5413 8251
Fax (81) 3 5413 8153

Germany
Tel (49) 257 1217 1200
Fax (49) 257 1207 1209

United States
Tel (877) 232-8095
Fax (956) 812-8086

553863-0001D TDS
Rev. 004 6/00

[7]

REFERENCES (Continued)

6. Maly, P., A.D. Thall, B. Petryniak, C.E. Rogers, P.L. Smith, R.M. Marks, R.J. Kelly, K.M. Gersten, G. Cheng, T.L. Saunders, S.A. Camper, R.T. Camphausen, F.X. Sullivan, Y. Isogai, O. Hindsgaul, U.H. von Andrian, and J.B. Lowe. 1996. The $\alpha(1,3)$ fucosyltransferase Fuc-TVII controls leukocyte trafficking through an essential role in L-, E-, and P-selectin ligand biosynthesis. *Cell* 86: 643 - 653.
7. Butcher, E.C. Personal communication.
8. Blinn, R.M., A. Whyte, S.T. Licence, A.A. Harrison, Y.T.M. Tsang, D.O. Haskard, and M.K. Robinson. 1996. The role of E-selectin in lymphocyte and polymorphonuclear cell recruitment into cutaneous delayed hypersensitivity reactions in sensitized pigs. *J. Immunol.* 157: 4094 - 4099.
9. DiCorlevo, T.G., K.D. Puri, R.A. Warnock, T.A. Springer, and U.H. von Andrian. 1996. Platelet-mediated lymphocyte delivery to high endothelial venules. *Science* 273: 252 - 255.
10. Puri, K.D., E.B. Finger, G. Gaudernack, and T.A. Springer. 1995. Sialomucin CD34 is the major L-selectin ligand in human tonsil high endothelial venules. *J. Cell Biol.* 131: 261 - 270.

For Research Use Only. Not For Diagnostic or Therapeutic Use.

Conditions: BD Pharmingen will not be responsible for violations or patent infringements which may occur with the use of our products.

Hazardous Ingredient: Sodium Azide. Avoid exposure to skin and eyes, ingestion, and contact with heat, acids, and metals. Wash exposed skin with soap and water. Flush eyes with water. Dilute with running water before discharge into plumbing.

BD Biosciences www.bdbiosciences.com

Biometric Imaging Cantech		Immunocytometry Systems Labware		Transduction Laboratories PharMingen	
Asia Pacific: Tel (852) 8610 423 Fax (852) 8601 390	Brazil Tel (55) 11 5185 0823 Fax (55) 11 5185 9834	Canada Tel (888) 258 0187 Fax (905) 542 8991 canada@bdb.com	Europe BD Biosciences Tel (32) 53 720 211 Fax (32) 53 720 430	Europe Life Science Research Tel (49) 6221 305 531 Fax (49) 6221 305 531	Japan Tel (81) 3 5413 8251 Fax (81) 3 5413 8155
				Madrid Tel (3291) 237 1200 Fax (3291) 237 1288	United States Tel 877 332 8995 Fax 858 612 8888

55380-088010-TDS
Rev 004 8/00

INTERNATIONAL SEARCH REPORT

International application No.

PCT/JP2004/009805

A. CLASSIFICATION OF SUBJECT MATTER

Int.Cl⁷ G01N33/574

According to International Patent Classification (IPC) or to both national classification and IPC

B. FIELDS SEARCHED

Minimum documentation searched (classification system followed by classification symbols)

Int.Cl⁷ G01N33/53-579

Documentation searched other than minimum documentation to the extent that such documents are included in the fields searched

Jitsuyo Shinan Koho 1922-1996 Toroku Jitsuyo Shinan Koho 1994-2004

Kokai Jitsuyo Shinan Koho 1971-2004 Jitsuyo Shinan Toroku Koho 1996-2004

Electronic data base consulted during the international search (name of data base and, where practicable, search terms used)

JICST (JOIS), CA (STN)

C. DOCUMENTS CONSIDERED TO BE RELEVANT

Category*	Citation of document, with indication, where appropriate, of the relevant passages	Relevant to claim No.
Y/A	UCHIMURA et al., "Specificities of N-acetyl glucosamine-6-0-sulfotransferases in relation to L-selectin ligand synthesis and tumor-associated enzyme expression" J.Biol.Chem., Vol.277, (2002), pages 3979 to 3984	7-9/1-6,10
Y/A	SEKO et al., "Ectopic expression of a GlcNAc 6-0-sulfotransferase, GlcNAc6ST-2, in colonic mucinous adenocarcinoma", Glycobiology, Vol. 12 (2002), pages 379 to 388	7-9/1-6,10
Y/A	JP 06-000093 A (Otsuka Pharmaceutical Co., Ltd.), 11 January, 1994 (11.01.94), (Family: none)	7-9/1-6,10

☒ Further documents are listed in the continuation of Box C.

☐ See patent family annex.

* Special categories of cited documents:

"A" document defining the general state of the art which is not considered to be of particular relevance

"E" earlier application or patent but published on or after the international filing date

"L" document which may throw doubts on priority claim(s) or which is cited to establish the publication date of another citation or other special reason (as specified)

"O" document referring to an oral disclosure, use, exhibition or other means

"P" document published prior to the international filing date but later than the priority date claimed

"T" later document published after the international filing date or priority date and not in conflict with the application but cited to understand the principle or theory underlying the invention

"X" document of particular relevance; the claimed invention cannot be considered novel or cannot be considered to involve an inventive step when the document is taken alone

"Y" document of particular relevance; the claimed invention cannot be considered to involve an inventive step when the document is combined with one or more other such documents, such combination being obvious to a person skilled in the art

"&" document member of the same patent family

Date of the actual completion of the international search
09 September, 2004 (09.09.04)

Date of mailing of the international search report
28 September, 2004 (28.09.04)

Name and mailing address of the ISA/
Japanese Patent Office

Authorized officer

Facsimile No.

Telephone No.

INTERNATIONAL SEARCH REPORT

International application No.

PCT/JP2004/009805

C (Continuation). DOCUMENTS CONSIDERED TO BE RELEVANT

Category*	Citation of document, with indication, where appropriate, of the relevant passages	Relevant to claim No.
Y/A	JP 06-046878 A (Tosoh Corp.), 22 February, 1994 (22.02.94), (Family: none)	7-9/1-6, 10
A	JP 11-313684 A (Seikagaku Corp.), 06 November, 1999 (06.11.99), & EP 943688 A & US 6037159 A	1-10

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